

Low pg/mL Detection of Underivatized 17 β -Estradiol in Serum Through Increased Ion Sampling Efficiency Using LC-MS/MS

Anabel Fandino¹, Brian C. Netzel², Ravinder J. Singh², ¹ Agilent Technologies Inc., ² Mayo Clinic, Rochester MN

ASMS 2011
WP 150



Introduction

17 β -estradiol (E2) is the subject of significant clinical research. Low levels of 17 β -estradiol have presented several challenges for traditional analysis of the molecule. Due to a lack of highly ionizable functional groups, mass spectrometric approaches have relied on laborious derivatization methods to achieve a sufficient lower limit of quantitation (LLOQ). Through the use of dual ion funnel technology, ion sampling efficiency has been improved to the point that underivatized 17 β -estradiol can be quickly and accurately quantified at low levels using an LC-MS/MS approach.

Experimental

Sample preparation

Calibration standards ranging from 1 to 500 pg/ml were prepared by dissolving E2 in charcoal stripped human serum. 500 μ L of spiked serum was extracted with diethyl ether. Samples were vortexed, centrifuged and 1 mL of the supernatant evaporated under nitrogen until dryness. Dried down extracts were reconstituted in water/methanol (70/30) and introduced to the LC-MS/MS system. In addition to these spiked serum samples, two sets of dried down extracts from the Mayo Clinic (Rochester, MN) were analyzed. The first set of samples consisted of phosphate buffer saline calibration standards in the range of 10 pg/mL to 600 pg/mL and patient samples. The second set consisted of calibration standards in the range of 0.1 pg/mL to 10 pg/mL.

Single LC Configuration: Agilent 1290 Infinity UHPLC


Column: Agilent Zorbax Extend-C18, 2.1 x 50 mm, 1.8 μ m. **Mobile phase:** A= 0.05% Ammonium Hydroxide in water, B= Acetonitrile; **Injection volume:** 10 μ L; **Flow rate:** 0.4 mL/min. **Gradient:** 15%B to 100%B in 1 min, 100%B during 1 min, 15%B at 2.01 min, stop time at 3.0 min. Autosampler with flexible cube

Matrix-Stripping LC Configuration: Agilent 1260 Infinity LC

Trapping Column: Agilent ZORBAX Extend-C18, 5 μ m, 4.6 x 12.5mm. **Analytical Column:** Agilent Zorbax Extend-C18, 2.1 x 50 mm, 1.8 μ m. **Injection V.:** 15 μ L or 100 μ L. **Loading:** 8%B. **Analytical Gradient:** 15%B to 100%B in 2 or 5 min, 100% B during 1 min. **Switching valve:** 0 min (Loading), 1 or 2 min. (Analysis), end of the gradient (Loading)

Experimental


Column-switching hardware and valve positions

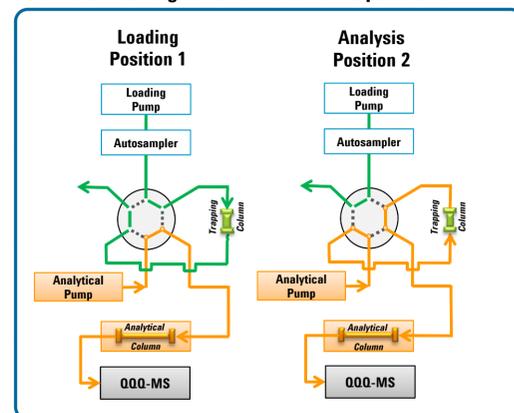


Agilent 6490 Triple Quadrupole MS:

Scan type: MRM. **Polarity:** Negative. **MRM transitions:** 271.3 \rightarrow 145.2 (Estradiol) and 276.3 \rightarrow 147.2 (Estradiol-d₅). **Collision energy:** 41 V. **Dwell time:** 100 ms. **Resolution:** MS1/MS2: Unit/Unit (Q1: 0.7 m/z / Q2: 0.7 m/z). **Agilent Jet Stream Parameters:** Drying gas temperature: 180 $^{\circ}$ C, Drying gas flow: 24 L/min, Sheath gas temperature: 350 $^{\circ}$ C, Sheath gas flow: 12 L/min, Nebulizer pressure: 25 psi, Nozzle voltage: 2000 V, Capillary voltage: 3000 V.


iFunnel Technology: Agilent Jet Stream (AJS) technology consists of the addition of a concentric super-heated nitrogen sheath gas to the nebulizer. The super-heated sheath gas collimates the nebulizer spray producing efficient desolvation and ion generation. The use of an hexabore capillary increases the interface area of the capillary inlet within the AJS thermal ion rich zone. The bores spread across the central, ion rich part of the AJS thermal confinement zone. In this way the multibore capillary captures more ions but also more gas. This gas load would normally overwhelm the vacuum system and therefore a novel dual ion funnel was added to efficiently remove the gas, while focusing the ions into the entrance of the first quadrupole.

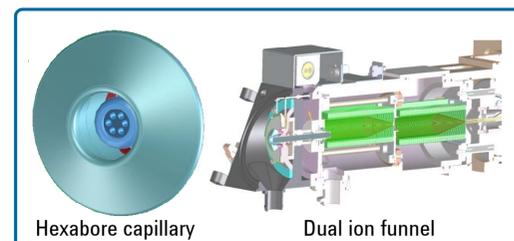


Single LC Configuration, Inj. V = 10 μ L: Assay performance

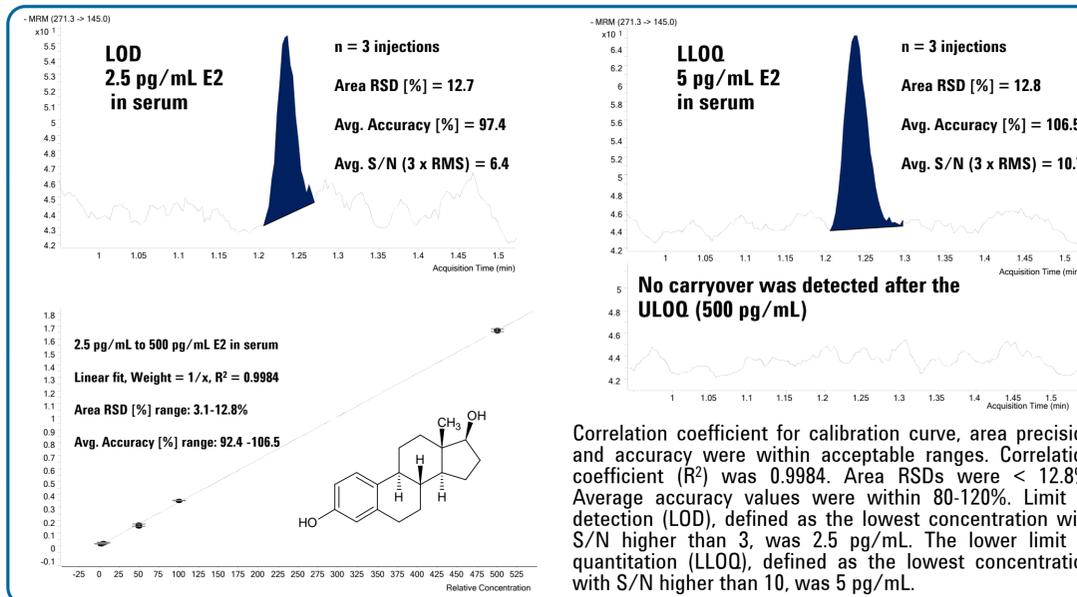

Matrix-Stripping LC Configuration, Inj. V = 15 μ L: Assay performance

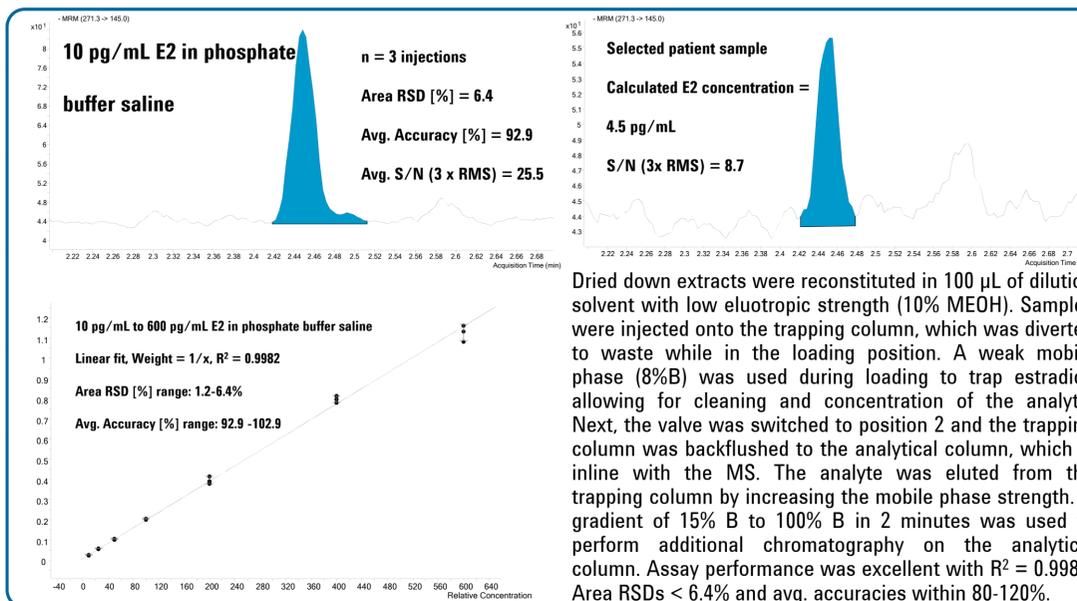

Matrix-Stripping LC Configuration, Inj. V = 100 μ L: Assay performance

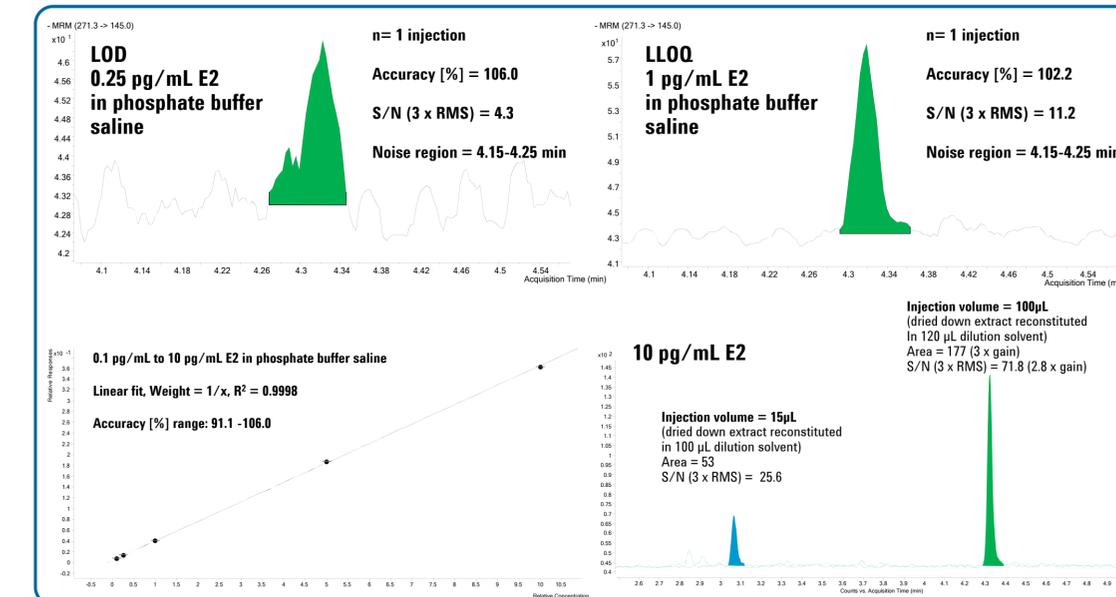

Although matrix-stripping LC involves more complicated column-switching hardware than is required for a single LC, it has several advantages. The setup allows for online sample cleanup, which reduces matrix effects, but it is also possible to inject large volumes of sample onto a narrow-bore sub-2 μ m particle size column without any band broadening effects. Therefore the possibility of achieving a LLOQ lower than 5 pg/mL was investigated by injecting 100 μ L onto the matrix-stripping LC/MS system. A theoretical area response gain factor of 5.5 was expected in comparison to a 15 μ L injection volume. Calibration standards in the range of 0.1 to 10 pg/mL E2 were prepared in phosphate buffer saline. Dried down extracts were reconstituted in 120 μ L dilution solvent (10% MEQH) and 100 μ L injected onto the LC/MS system. A longer loading time (2 min.) and longer gradient time (15%B to 100% in 5 min.) were used for the analysis of this set of samples. A factor 3.3 area response gain and factor 2.8 S/N (3 x RMS) gain were observed in comparison to the method using injection volume = 15 μ L. Correlation coefficient for calibration curve and accuracy were within acceptable ranges. Correlation coefficient (R²) was 0.9998. Accuracy values were within 80-120%. The LOD was 0.25 pg/mL and the LLOQ was 1 pg/mL.

Conclusions

- Taking full advantage of dual ion funnel technology, it was possible to quantify underivatized estradiol at the low pg/mL level
- Sensitivity of the method was demonstrated with LLOQ at 1 pg/mL underivatized E2 using a matrix-stripping LC approach capable of concentrating a 100 μ L injection volume to avoid band broadening effects
- This high quality estradiol assay with LLOQ = 1 pg/mL and excellent performance will be an important tool for clinical research where assay sensitivity is one of the most important factors