

Metabolic Stability Study Using Cassette Analysis and Polarity Switching in an Ultra High Performance Liquid Chromatography (UHPLC)-Triple Quadrupole LC/MS System

Application Note

Authors

Anabel S. Fandiño, Edgar Naegele
Agilent Technologies, Inc.
Santa Clara, CA USA

Abstract

This application note demonstrates the performance of the Agilent 1290 Infinity LC System coupled to an Agilent 6460A Triple Quadrupole LC/MS for the evaluation of metabolic stability using the cassette analysis approach combined with fast polarity switching.

We present a UHPLC/MS/MS method that demonstrates:

- High throughput analysis (retention times less than 0.8 min at a mobile phase flow rate of 1.5 mL/min and binary pump pressure 1,100 bar)
- High speed MS/MS analysis with cycle time as low as 114 ms to monitor 6 transitions and switch polarity (polarity switching time 30 ms).
- Good precision with area RSD [%] or relative area RSD [%] less than 10%, enabled by collecting enough data points (at least 9 points) across the chromatographic peaks (as narrow as 0.37 sec at half height).
- Good accuracy over the expected concentration range using either external or internal standard calibration (within 90.7 – 107.1%)
- Comparable results (% parent drug remaining) when using polarity switching or non-switched analysis
- Excellent chromatographic resolution to avoid overestimation in quantitation due to in-source CID conversion of a metabolite into its parent drug



Agilent Technologies

Introduction

In vitro ADME assays like metabolic stability play an important role in the early understanding of *in vivo* pharmacokinetic characteristics and help to discard non-drug-like compounds that would fail later stages of development. Thus, fast screening methods are needed to screen the large collection of new chemical entities that need characterization in the early stages of drug discovery. Typically, metabolic stability assessment is performed using the cassette approach and fast LC/MS/MS methods. In the cassette approach, a cocktail of substrates ionizing either in positive or in negative mode is created post-incubation in order to reduce the number of samples to analyze. Fast LC methods, which use high mobile phase flow rates and short sub-2 μm columns, allow reducing the analysis time and therefore increasing sample throughput [1, 2]. Detection is usually performed using a triple quadrupole mass spectrometer working in MRM scan mode.

Cassettes are typically designed so that substrates ionizing either in positive or negative mode are pooled together for LC/MS/MS analysis. The reason for that is because MS/MS analysis may show challenges with quantitative data quality if the MS detector is not able to acquire data fast enough in order to collect at least 9 to 10 data points across the extremely narrow peaks generated using fast chromatography. Nevertheless, pooling of substrates ionizing in positive and negative mode would present a pronounced flexibility advantage.

This work describes the advantage of combining not only cassette analysis and fast LC but also using fast polarity switching MS/MS analysis to increase throughput and flexibility in a metabolic stability assay, while maintaining good precision and accuracy. Moreover, the Agilent 1290 Infinity LC system flow range of up to 2 mL/min at 1,200 bar allowed excellent chromatographic resolution between critical peaks.

Experimental

Reagents and supplies

Buspirone hydrochloride and verapamil hydrochloride were purchased from Sigma-Aldrich (St. Louis, MO USA). Dextromethorphan, diclofenac, and dextrorphan- d_3 (ISTD) were purchased from Cerilliant (Round Rock, TX USA). Diclofenac- d_4 (phenyl- d_4) was purchased from CDN Isotopes (Quebec, Canada). NADPH regeneration system, solution A (beta-nicotinamide adenine dinucleotide phosphate sodium salt (NADP), glucose-6-phosphate and MgCl_2 in water) and solution B (glucose-6-phosphate dehydrogenase in 5 mM sodium citrate) as well as 0.5 M potassium phosphate buffer, pH 7.4 were purchased from BD Biosciences (Woburn, MA USA). Male Sprague Dawley rat liver S9 homogenate was purchased from In vitro Technologies (Baltimore, MD). All other reagents and organic solvents were of analytical grade and from VWR (Darmstadt, Germany).

Incubation of substrates with rat liver S9 fractions

The incubation mixtures for phase I metabolism consisted of an amount of S9 preparation equivalent to 0.3 mg protein, 1 μM substrate (buspirone, verapamil, dextromethorphan or diclofenac from a 100 μM stock solution in 70% water / 30% acetonitrile), 1.3 mM NADP, 3.3 mM glucose-6-phosphate, 3.3 mM magnesium chloride and 0.4 U/mL glucose-6-phosphate dehydrogenase in 0.1 M phosphate buffer

(pH 7.4) made up to a total volume of 300 μL . Incubation was carried out at 37°C. A cocktail of the 4 substrates was created post-incubation to reduce analysis time. A 25 μL aliquot was taken at 0, 5, 10, 15, 25 and 35 min from each incubate and the reaction was stopped by adding 300 μL acetonitrile containing the internal standards (dextrorphan- d_3 and diclofenac- d_4) followed by centrifugation for 10 min at 14,000 g. The supernatant was evaporated to dryness using a gentle stream of nitrogen and reconstituted with water/acetonitrile 80/20 v/v containing 0.1% formic acid for UHPLC/MS/MS analysis.

Equipment

- Agilent 1290 Infinity LC System comprising 1290 Infinity Binary Pump with integrated degasser, 1290 High Performance Autosampler with Thermostat and 1290 Infinity Thermostatted Column Compartment
- Agilent 6460A Triple Quadrupole LC/MS System with Agilent Jet Stream technology
- Agilent MassHunter Workstation software for instrument control, data acquisition and data processing
- Agilent MassHunter Optimizer Software
- Agilent Rapid Resolution High Definition (RRHD) Zorbax SB-C18, 2.1 x 50 mm, 1.8 μm column

Agilent 1290 Infinity Methods

The mobile phase consisted of:

Solvent A: water with 0.1% formic acid

Solvent B: acetonitrile with 0.1% formic acid

Injection volume: 1 µL

Needle wash: 20 sec in flushport with methanol/water 50/50 v/v (0.1% formic acid)

UHPLC Method 1:

Column temperature: 25°C or 40°C

Flow rate: 1.0 mL/min

Gradient: 25% B during 0.2 min, 80% B at 1 min, 80% B at 1.25 min, 25% B at 1.26 min, stop time at 1.8 min.

UHPLC Method 2:

Column temperature: 60 °C

Flow rate: 1.5 mL/min

Gradient: 25% B during 0.2 min, 80% B at 0.73 min, 80% B at 1.00 min, 25% B at 1.01 min, stop time at 1.5 min.

Agilent 6460A triple quadrupole conditions

Scan type: MRM (MassHunter optimizer software allows to quickly determine the optimal fragmentor voltage, MRM transitions and collision energy for the selected transitions. Optimized parameters are shown in **Tables 1, 2 and 3**)

Polarity: positive/negative, positive only or negative only

Parameters: drying gas temperature: 350°C, drying gas flow: 10 L/min, sheath gas temperature: 400°C, sheath gas flow: 12 L/min, nebulizer pressure: 35 psig, nozzle voltage: 0 V (+) 1000 V (-), capillary voltage: 4,000 V (+/-)

Polarity switching time: 30 ms

Compound	ISTD	Precursor Ion	MS1 Res	Product Ion	MS2 Res	Dwell time [ms]	Frag V [V]	CE [V]	Polarity
Verapamil		455.3	Unit	165.1	Unit	5	185	26	Positive
Buspirone		386	Unit	122.1	Unit	5	185	28	Positive
Diclofenac -d ₄	X	298	Unit	254	Unit	5	80	4	Negative
Diclofenac		294	Unit	250	Unit	5	95	5	Negative
Dextromethorphan		272.2	Unit	215.2	Unit	5	190	22	Positive
Dextrorphan -d ₃	X	261.2	Unit	157.1	Unit	5	125	40	Positive

Table 1: MRM acquisition with positive/negative switching.

Compound	ISTD	Precursor Ion	MS1 Res	Product Ion	MS2 Res	Dwell time [ms]	Frag V [V]	CE [V]	Polarity
Verapamil		455.3	Unit	165.1	Unit	22	185	26	Positive
Buspirone		386	Unit	122.1	Unit	22	185	28	Positive
Diclofenac		296	Unit	250	Unit	22	80	8	Positive
Dextromethorphan		272.2	Unit	215.2	Unit	22	190	22	Positive
Dextrorphan -d ₃	X	261.2	Unit	157.1	Unit	22	125	40	Positive

Table 2: MRM acquisition in positive polarity.

Compound	ISTD	Precursor Ion	MS1 Res	Product Ion	MS2 Res	Dwell time [ms]	Frag V [V]	CE [V]	Polarity
Arbitrary transition 1		455.3	Unit	165.1	Unit	22	185	26	Negative
Arbitrary transition 2		386	Unit	122.1	Unit	22	185	28	Negative
Diclofenac		294	Unit	250	Unit	22	95	5	Negative
Arbitrary transition 3		272.2	Unit	215.2	Unit	22	190	22	Negative
Diclofenac -d ₄	X	298	Unit	254	Unit	22	80	4	Negative

Table 3: MRM acquisition in negative polarity.

Results and Discussion

Speed and data quality using fast polarity switching

Operating the UHPLC system at high flow rates (1.0 mL/min or 1.5 mL/min) and pressures up to 1,100 bar enabled a run time of only 1.5 min and generated peak widths less than a second at half height (typically 0.4 to 1.0 sec). Sufficient data points across the peak (> 9 points) could be collected due to

low MS cycle times, ensuring high precision quantitation (area RSD [%] and relative area RSD [%] = area target / area internal standard < 10). **Fig. 1** shows in comparison the MRM chromatograms, data quality and analysis time achieved using flow rates of 1.5 and 1.0 mL/min.

Good quantitation requires 9 – 10 data points across a peak. As peaks get narrower, the MS detector must be able to acquire faster.

MS cycle times are reduced. During each cycle, the MS system must monitor 6 transitions and switch polarity. This is enabled by 30 ms polarity switching time and a total cycle time of only 114 ms as shown in **Fig. 2** for the substrate diclofenac. The graphic shows that even for peaks exhibiting a width of 0.37 sec at half height, the mass spectrometer was able to collect a total of 9 points, ensuring good area and relative area precision of 4.8 and 7.9%, respectively.

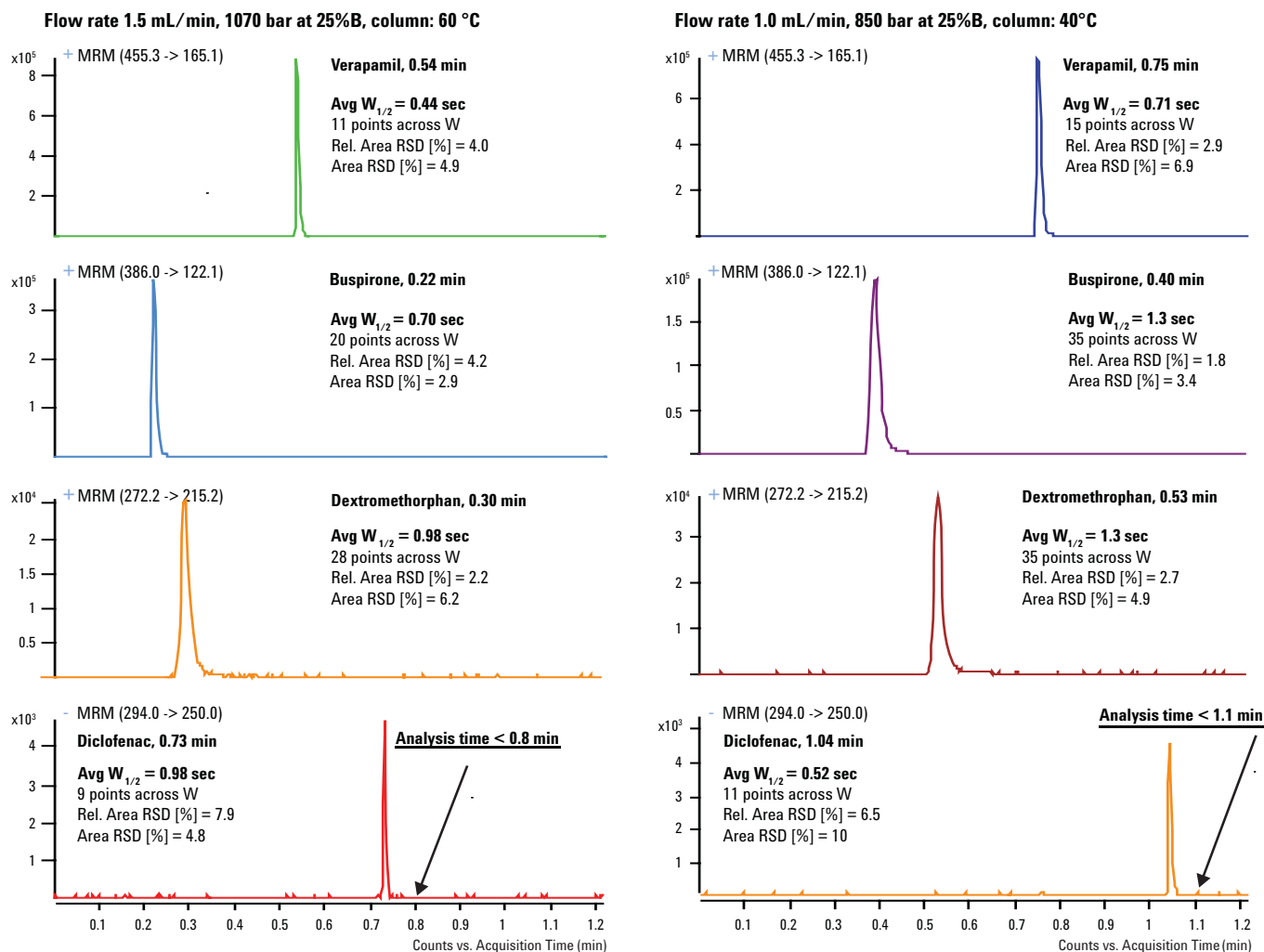


Fig. 1: MRM chromatograms, data quality and analysis time achieved using flow rates of 1.5 and 1.0 mL/min.

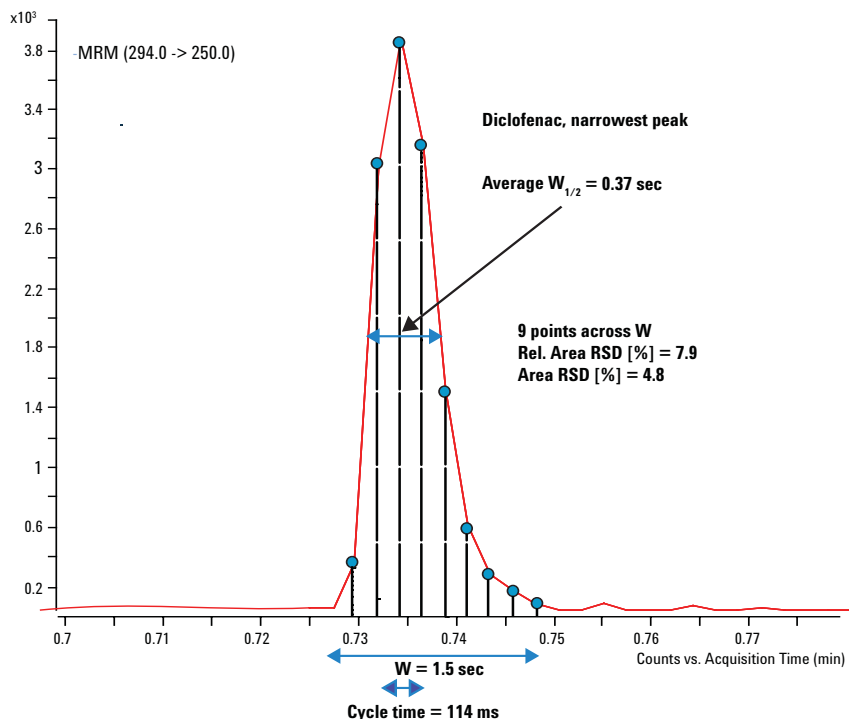


Fig. 2: MRM chromatogram of diclofenac obtained using a flow rate of 1.5 mL/min. The graphic shows the peak width at half height (0.37 sec) and the number of data points collected across the peak (9 data points).

Evaluation of linearity, precision and accuracy using fast polarity switching

Linearity, precision and accuracy were tested using standards at a flow rate of 1.0 mL/min and column temperature 25°C (UHPLC method 1). The % of remaining parent was simulated over the expected range. The results obtained from the linearity test using either external or internal calibration demonstrate that the cassette analysis with polarity switching approach can provide linear results over the desired range. Moreover, good precision and accuracy values were achieved at all levels as shown in **Table 4**. When using internal standard calibration, a conflict between $^{35}\text{Cl}_2$ -Diclofenac- d_4 (m/z 298) and $^{37}\text{Cl}_2$ -Diclofenac (m/z 298) was observed. The response of the transition 298 → 254 arising from $^{37}\text{Cl}_2$ -Diclofenac contributes to the response of the internal standard $^{35}\text{Cl}_2$ -Diclofenac- d_4 . This led to a quadratic curve due to artificially “increasing” internal standard. To solve this, either external calibration or a unique transition to diclofenac- d_4 like 298 → 217 (CE = 17) can be used.

	Buspirone		Verapamil		Dextromethorphan		Diclofenac	
	external calibration	isotope dilution	external calibration	isotope dilution	external calibration	isotope dilution	external calibration	isotope dilution
Type	Linear	Linear	Linear	Linear	Linear	Linear	Linear	Quadratic
r^2	0.9989	0.9994	0.9965	0.9989	0.9978	0.9980	0.9990	0.9991
Accuracy range [%]	97.4 – 102.0	94.4 – 107.1	90.7 – 105.1	96.2 – 103.2	97.6 – 102.4	95.5 – 106.1	96.2 – 103.1	99.9 – 100.4
Area RSD [%]	2.4 – 6.1		2.3 – 6.9		3.7 – 5.1		1.7 – 8.0	
Rel. area RSD range [%]	n.a.	1.7 – 4.9	n.a.	1.4 – 6.0	n.a.	2.6 – 5.9	n.a.	0.9 – 4.1

Table 4: Calibration curve type, Correlation coefficient (r^2), accuracy range [%], area RSD [%] and relative area RSD range [%] obtained using cassette analysis and fast polarity switching. References: n.a.: not applicable


Switched versus non-switched analysis – Area response and precision for pooled incubates

Average area responses obtained using switched analysis (fast polarity switching) or non-switched analysis (positive only or

negative only analysis) were similar. Area RSD and relative area RSD values were less than 10.4% using polarity switching and less than 8.6% using non-switched analysis.

Fig. 3 shows in comparison the MRM

chromatograms, average area response and relative area RSD [%] obtained for the pooled incubates analyzed using a flow rate of 1.0 mL/min and column temperature 25°C (UHPLC method 1).

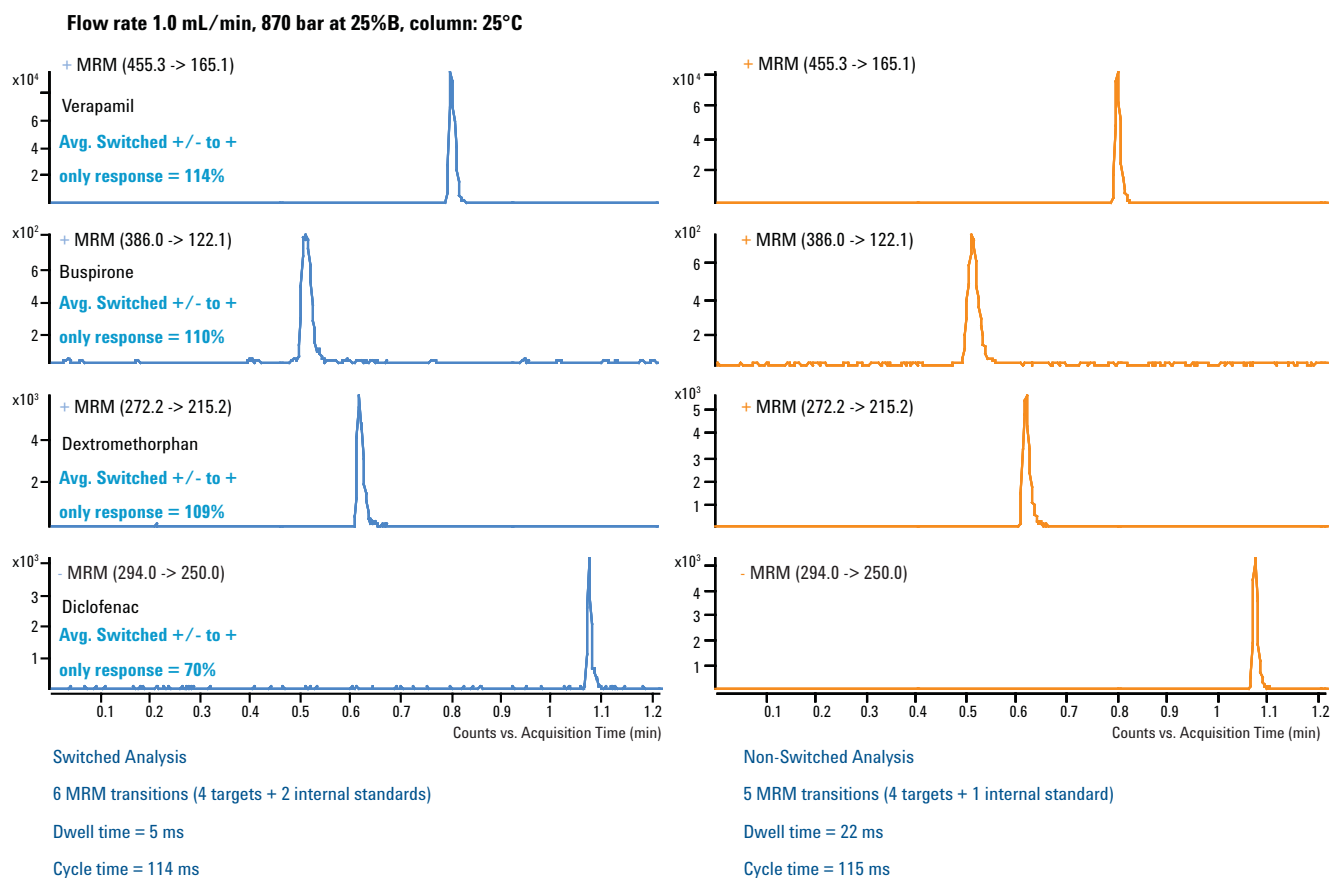

Time [min]	Buspirone		Verapamil		Dextromethorphan		Diclofenac			
	pos./neg.	pos. only	pos./neg.	pos. only	pos./neg.	pos. only	pos./neg.	neg. only	pos./neg.	neg. only
0	2.7	0.5	8.7	1.0	3.0	2.7	3.2	1.5	5.3	3.6
5	5.2	3.3	8.1	2.8	5.7	2.9	6.4	1.4	3.9	7.0
10	7.2	4.9	6.5	5.5	5.0	3.7	8.1	2.5	2.4	4.1
15	4.0	3.9	4.8	2.0	4.1	2.4	8.9	3.0	10.2	5.2
25	3.2	4.3	5.6	1.7	10.4	1.7	6.8	2.4	4.3	4.8
35	9.2	3.8	5.4	5.8	8.9	8.6	6.3	2.9	6.9	8.6

Fig. 3: MRM chromatograms showing average area response obtained using fast polarity switching (on the left) in comparison to non-switched analysis (on the right). The table at the bottom shows the precision obtained using fast polarity switching or non-switched analysis for all incubation times.

Switched versus non-switched analysis – Metabolic stability (% parent drug remaining)

The % parent drug remaining was determined by comparing the average relative area of the parent compound measured in

the 5-, 10-, 15-, 25- and 35-min samples to that in the 0-min sample and was calculated as follows:

$$\% \text{ parent remaining} = (\text{avg. rel. area at } t_x / \text{avg. rel. area at } t_0) \times 100$$

Fig. 4 shows the correlation of results obtained for verapamil using fast polarity switching and non-switched analysis. The results of the switched and non-switched analysis were comparable for all substrates with $r^2 > 0.9848$.

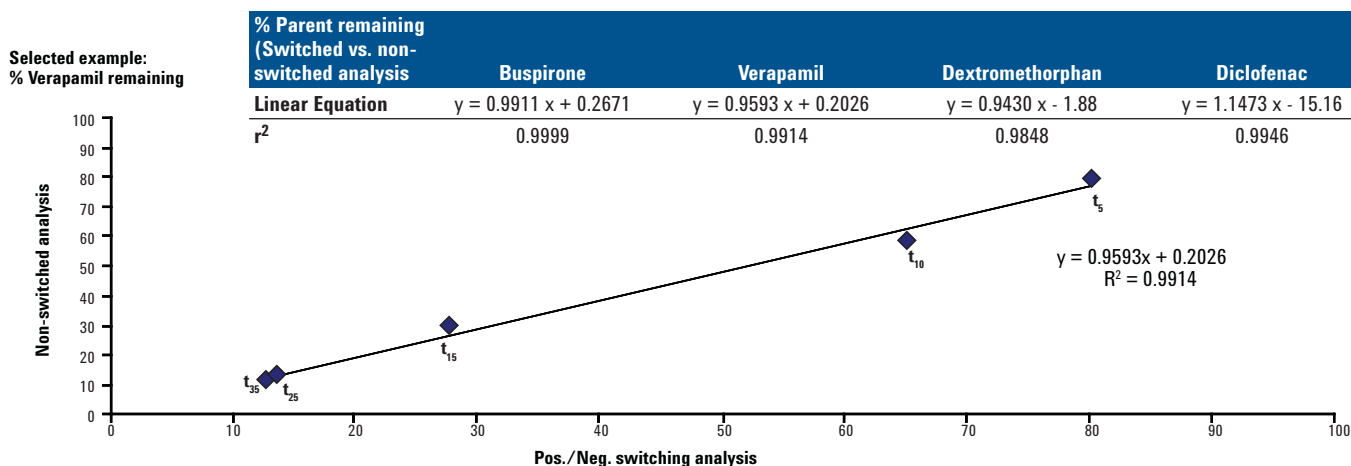


Fig. 4: Comparison of % parent drug remaining from fast polarity switching with non-switched analysis using UHPLC method 1 (flow rate = 1.0 mL/min, column temperature = 25°C).

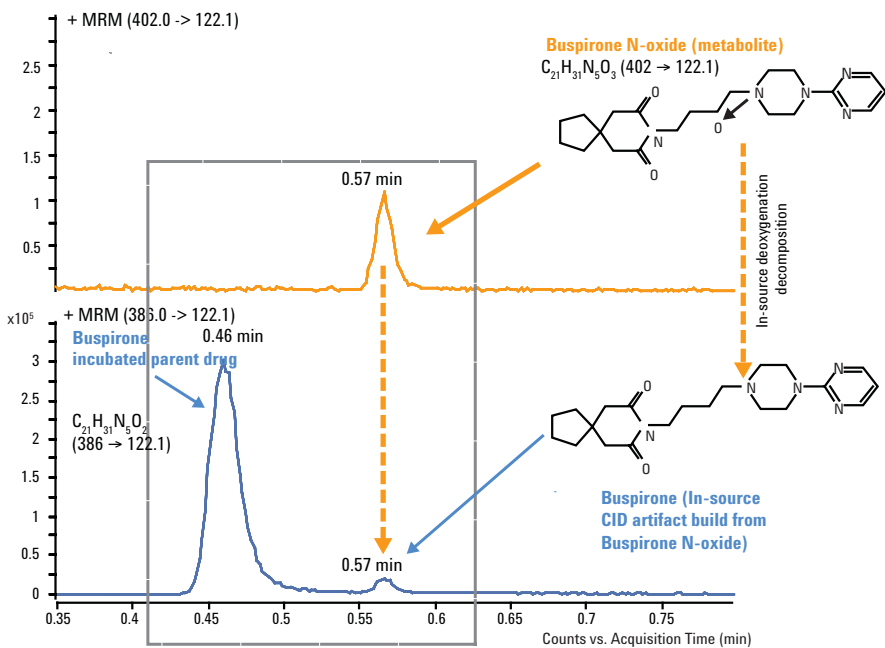


Fig. 5: Excellent chromatographic resolution achieved between the incubated parent drug buspirone and one of its metabolites (buspirone N-oxide) which builds buspirone by in-source deoxygenation.

Chromatographic resolution – Advantage for cassette analysis

Excellent resolution between critical peaks was demonstrated. During the determination of metabolic stability, it is important to achieve high chromatographic resolution to ensure that the metabolites produced by the hepatic S9 fraction are resolved from the parent compounds. Coelution may result in overestimation of the amount of parent compound present if the metabolite is thermally unstable and can convert into the parent compound by in-source CID; **Fig. 5** shows as an example the conversion of buspirone N-oxide to buspirone (parent drug) by deoxygenation. Excellent resolution between buspirone and buspirone N-oxide allowed to distinguish between the parent drug buspirone (originally present in the incubation sample, retention time = 0.46 min) and the artifactually built buspirone (built by in-source CID from buspirone N-oxide, retention time = 0.56 min).

Conclusions

The combination of cassette analysis and fast polarity switching in the Agilent UHPLC/MS/MS system provided:

- Flexibility for cassette design: Compounds ionizing in positive and negative mode can be pooled together.
- High throughput analysis (retention times less than 0.8 min at a mobile phase flow rate of 1.5 mL/min and binary pump pressure of 1,100 bar)
- High speed MS/MS analysis with cycle time as low as 114 ms to monitor 6 transitions and switch polarity (polarity switching time 30 ms).
- Good precision with area RSD [%] and relative area RSD [%] less than 10%, enabled by collecting enough data points (at least 9 points) across the chromatographic peaks (as narrow as 0.37 sec at half height).
- Good accuracy over the expected concentration range using either external or internal standard calibration (accuracy within 90.7 - 107.1%)
- Comparable results (% parent drug remaining) when using polarity switching or non-switched analysis
- Excellent chromatographic resolution to avoid overestimation in quantitation due to in-source CID conversion of a metabolite into its parent drug

References

1. Timothy J. Carlson and Michael B. Fisher, "Recent Advances in High Throughput Screening for ADME Properties," *Combinatorial Chemistry and High Throughput Screening*, Vol 11, No 3, March 2008, 258-264.
2. Walter A. Korfmacher, "Principles and applications of LC-MS in drug discovery", *Drug Discovery Today*, Vol 10, No 20, October 2005, 1357-1367.

www.agilent.com/chem/qqq

This item is intended for Research Use Only. Not for use in diagnostic procedures. Information, descriptions, and specifications in this publication are subject to change without notice.

Agilent Technologies shall not be liable for errors contained herein or for incidental or consequential damages in connection with the furnishing, performance or use of this material.

© Agilent Technologies, Inc. 2009

Published in the U.S.A., September 15, 2009
5990-4469EN



Agilent Technologies