



# C1 – C3 Hydrocarbon Analysis Using the Agilent 490 Micro GC – Separation Characteristics for PoraPLOT U and PoraPLOT Q Column Channels

## Application Note

Micro Gas Chromatography, Hydrocarbon analysis

### Author

Remko van Loon  
Agilent Technologies, Inc.  
Middelburg  
The Netherlands



### Introduction

This application note shows the possibilities and limitations in fast analysis of saturated and unsaturated C1 to C3 hydrocarbons using an Agilent 490 Micro GC. The chromatograms and results outline the similarities and differences when using a PoraPLOT U and a PoraPLOT Q column channels. Both the PoraPLOT U and the PoraPLOT Q are capable of resolving methane from the composite air peak and separate CO<sub>2</sub> from methane and the C2 hydrocarbons.

The PoraPLOT U column channel will have the following separation characteristics:

- Baseline separation for ethane, ethylene and acetylene
- Coelution of propane and propylene

The separation characteristics for the PoraPLOT Q column channel are:

- Coelution of ethylene and acetylene
- Baseline separation for propane and propylene





**Agilent Technologies**

If you want to the ability to measure anywhere and get the results you need in seconds, the Agilent 490 Micro GC is the ideal solution. With its rugged, compact, laboratory quality gas analysis platform, the 490 Micro GC generates more data in less time for faster, and better, business decisions.

## Instrumentation

For this application an Agilent 490 Micro GC (G3581A) equipped with a PoraPLOT U and a PoraPLOT Q was used. The setup parameters for the column is found in the table below.

	PoraPLOT U, 10 m	PoraPLOT Q, 10 m
Column temperature	80 °C	80 °C
Carrier gas	Helium, 200 kpa	Helium, 200 kpa
Injector temperature	110 °C	110 °C
Injection time	20 ms	20 ms

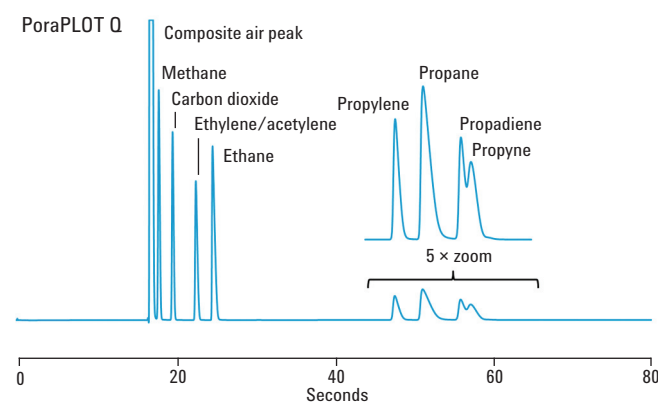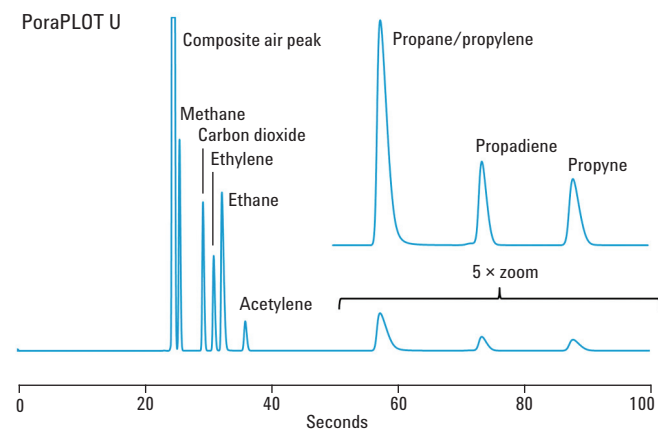

## Sample information

Nitrogen	Balance
Methane	5.0 %
Carbon dioxide	3.0 %
Ethylene	2.0 %
Ethane	4.0 %
Acetylene	1.0 %
Propylene	1.0 %
Propane	2.0 %
1,2-Propadiene	0.97 %
Propyne	0.99 %

## For More Information

These data represent typical results. For more information on our products and services, visit our Web site at [www.agilent.com/chem](http://www.agilent.com/chem).

[www.agilent.com/chem](http://www.agilent.com/chem)

Agilent shall not be liable for errors contained herein or for incidental or consequential damages in connection with the furnishing, performance, or use of this material.

Information, descriptions, and specifications in this publication are subject to change without notice.

© Agilent Technologies, Inc., 2012  
 Printed in the USA  
 January 11, 2012  
 5990-9165EN



**Agilent Technologies**