

Reveleris® Flash Cartridges





Get the most out of your flash applications by choosing Reveleris® flash cartridges, the only flash cartridges designed by a silica expert for ideal chromatographic performance.

- Increased loading capacity
- Durable 200psi design
- Superior resolution with ultra-pure 40µm silica
- Four unique chemistries available



Luer design is compatible with most instruments, and adaptors are available for non-luer connections.


Flow Rate Comparison in the Purification of α from δ -tocopherol

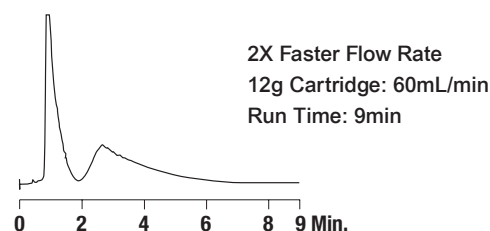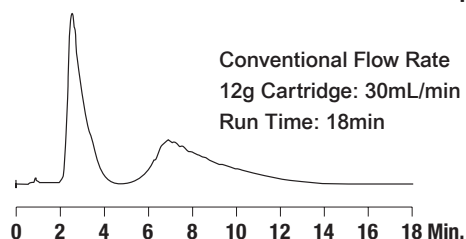

Flash Cartridge: Reveleris® 12 gram cartridge
Wavelength: 295nm
Isocratic: Solvent A: Hexane 97%
Solvent B: THF 3%
Sample: 500mg/mL
Sample Inject Volume: 2mL

Higher resolution of Reveleris® cartridges enables the flow rate to be increased 2X while maintaining base-line separation. For this application, that results in over a 40% time savings per run.

Flavanone Glycoside Purification

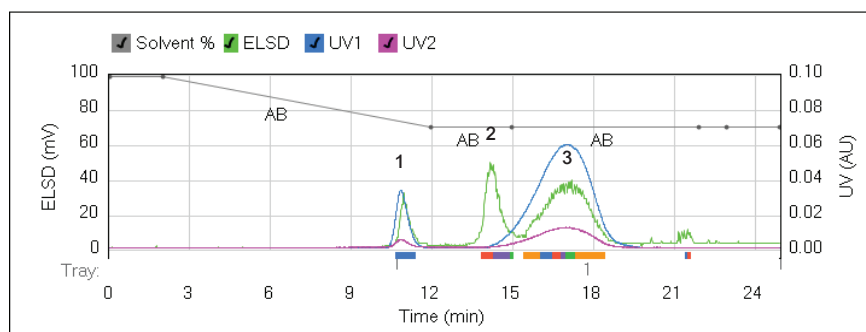

Column: Reveleris® Amino 40g flash cartridges
Flow Rate: 40mL/min
Equilibration: 10.0min
Run Length: 25.0min

UV1 Wavelength: 220nm
UV2 Wavelength: 254nm
ELSD Carrier: Iso-propanol
Solvent A: Water
Solvent B: Acetonitrile

1. Hesperitin 20mg
2. Glucose 20mg
3. Neohesperidin 20mg

Gradient Table

Step	Min	Solvents	% 2nd
1	0.0	AB	99
2	2.0	AB	99
3	10.0	AB	70
4	13.0	AB	70

Combining the high capacity of the amino phase cartridge and RevealX™ detectoin technology, neohesperidin is isolated and purified from glucose.

Visit www.discoverysciences.com/FlashChromatography/ApplicationNotes.aspx to browse the latest applications utilizing Reveleris® flash cartridges



Unique Chemistries for Challenging Applications

While silica is the most popular choice for flash applications, other chemistries bring unique properties that can help resolve complex samples.

Chemistry	Ideal For
Silica	Conventional applications, polar compounds
Reversed Phase (C18)	Non-polar compounds
Amino	Better peak shapes with polar and basic compounds; ideal for carbohydrates
Diol	Alternate selectivity for a variety of polarity types; ideal for lipids

Reveleris® Flash Cartridges

Description	Qty	Part No.
Silica		
4g Silica Cartridge	20	5146130
12g Silica Cartridge	20	5146131
40g Silica Cartridge	15	5146132
80g Silica Cartridge	12	5146133
120g Silica Cartridge	10	5146134
330g Silica Cartridge	4	5146135
Reversed Phase (C18)		
4g C18 Cartridge	each	5152990
12g C18 Cartridge	each	5152103
40g C18 Cartridge	each	5152104
80g C18 Cartridge	each	5152105
120g C18 Cartridge	each	5152991
Amino		
4g Amino Cartridge	each	5157330
12g Amino Cartridge	each	5157331
40g Amino Cartridge	each	5157332
80g Amino Cartridge	each	5157333
120g Amino Cartridge	each	5157334
Diol		
4g Diol Cartridge	each	5157335
12g Diol Cartridge	each	5157336
40g Diol Cartridge	each	5157337
80g Diol Cartridge	each	5157338
120g Diol Cartridge	each	5157339

Application Kits for Simplified Method Development

Reveleris® Application Kits combine a variety of cartridges tailored to different application areas. Rapidly select the best chemistry for many applications and compound classes.

Application Kits

Description	Part No.
Radiochem Flash Kit	5157384
Carbohydrate Flash Kit	5157385
Impurity Isolation Flash Kit	5157386
Lipid Flash Kit	5157387
Natural Products Flash Kit	5157388

Screening Columns for Efficient Planning

Reveleris® C18 HPLC columns are ideal for prescreening by HPLC prior to flash purification.

Reveleris® C18 HPLC Columns*

Description	Qty	Part No.
LC Column C18, 4.6 x 50mm	3	5156381
LC Column C18, 4.6 x 150mm	3	5156382
LC Column C18, 4.6 x 250mm	3	5156383

*For use with Reveleris® Navigator Feature

Accessories

Description	Qty	Part No.
Cartridge Adaptor		
Flash Cartridge Adaptors Kit		5156387
<i>Use with 4-330g cartridges</i>		
Reveleris® TLC Plates		
Silica 20x20cm	25	8621250
C18 10x20	25	8625320
Adjustable Solid Loader Plunger		
3g (15mL)	each	5148510
15g (75mL)	each	5148511
40g (150mL)	each	5148512
Empty Solid Loaders		
15mL	100	5142050
75mL	100	5142051
150mL	50	5147808
Filled Solid Loaders		
3g (15mL)	20	5152371
15g (75mL)	16	5142034
40g (150mL)	12	5147772



For maximum performance, use Reveleris® flash cartridges with the Reveleris® flash instrument. Take advantage of higher pressures, user-friendly software, advanced detection capabilities, and RFID tagging. The Navigator software feature automatically converts HPLC and TLC prescreening data into optimized flash conditions. Contact us or visit our website to learn more.

www.discoverysciences.com

ALLTECH® and REVELERIS® are trademarks, registered in the United States and/or other countries, of Alltech Associates, Inc. REVEALX™ is a trademark of Alltech Associates, Inc. GRACE®, and GRACE DAVISON® are trademarks, registered in the United States and/or other countries, of W. R. Grace & Co.-Conn. GRACE DAVISON DISCOVERY SCIENCES™ and GRACERESOLV™ are trademarks of W. R. Grace & Co.-Conn. This trademark list has been compiled using available published information as of the publication date of this brochure and may not accurately reflect current trademark ownership or status. Alltech Associates, Inc. is a wholly owned subsidiary of W. R. Grace & Co.-Conn. Grace Davison Discovery Sciences is a product group of W. R. Grace & Co.-Conn., which now includes all product lines formerly sold under the Alltech brand. © Copyright 2011 Alltech Associates, Inc. All rights reserved. The information presented herein is derived from our testing and experience. It is offered for your consideration and verification. Since operating conditions vary significantly, and are not under our control, we disclaim all warranties on the results that may be obtained from the use of our products. W. R. Grace & Co.-Conn. and its subsidiaries can not be held responsible for any damage or injury occurring as a result of improper installation or use of its products. Grace reserves the right to change prices and/or specifications without prior notification.

01/2011 M372

W. R. Grace & Co.-Conn. 7500 Grace Drive, Columbia, MD 21044

Grace Davison Discovery Sciences Regional Headquarters:

IN THE AMERICAS:
2051 Waukegan Road
Deerfield, IL 60015-1899
Tel: 1 847 948 8600
Fax: 1 847 948 1078
Email: discoverysciences@grace.com

IN EUROPE:
Brandstraat 12
B-9160 Lokeren, Belgium
Tel: +32 (0)9 340 65 65
Fax: +32 (0)9 340 65 60
Email: discoverysciences.EU@grace.com

IN ASIA:
19th Floor, K.Wah Center
1010 Huai Hai Zhong Road
Shanghai, 200031 PRC
Tel: 86 21 5467 4678
Fax: 86 21 5405 1500
Email: dsbiz.asia@grace.com

IN INDIA:
17 Commerce Center
Opposite Krishna Hospital, Paud Road
Pune, 411038 India
Tel: +91 20 6644 9900
Fax: +91 20 2544 1740
Email: pune@grace.com

IN AUSTRALIA/NEW ZEALAND:
2 Kerr Court
Rowville, 3178
Victoria, Australia
Tel: +61 3 9237 6100
Fax: +61 3 9237 6101
Email: discoverysciences.AU@grace.com

GRACE